

Mouse IgA ELISpot^{BASIC}

Product Code: 3835-2A

CONTENTS:

Vial 1 (blue top)

Anti-IgA capture antibody (800 μ l)

Concentration: 0.5 mg/ml

Vial 2 (red top)

Anti-IgA biotinylated detection antibody (100 μ l)

Concentration: 0.5 mg/ml

Vial 3 (white top)

Streptavidin-Alkaline Phosphatase (50 μ l)

Vial 4 (black top)

Polyclonal activator: R848 (100 μ l)

Concentration: 1 mg/ml

Vial 5 (blue top)

Lyophilised recombinant mouse IL-2 (0.5 μ g)


Antibodies are supplied in sterile filtered (0.2 μ m) PBS with 0.02% sodium azide. Streptavidin-ALP is supplied in 0.1 M Tris buffer with 0.002% Kathon CG. R848 is supplied in sterile filtered (0.2 μ m) PBS containing 1% DMSO. Vials have been overfilled to ensure recovery of the specified amount.

STORAGE:

Shipped at ambient temperature. On arrival antibodies and Streptavidin-ALP should be stored refrigerated at 4-8°C. R848 and IL-2 should be stored at -20°C.


Schematic illustration of the different protocols

PROTOCOL I

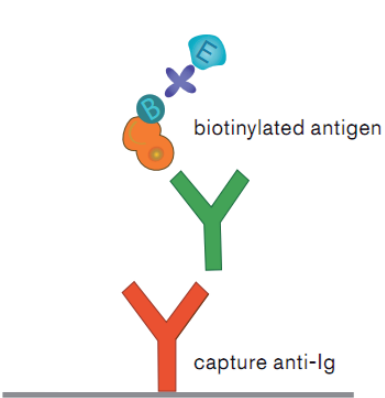

Antigen-specific Ig ELISpot

PROTOCOL II

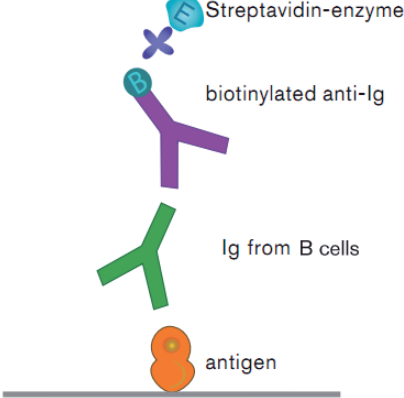

Total Ig ELISpot

Kit description

Please read through before starting the assay

This kit contains reagents for detection and enumeration of B cells secreting IgA antibodies. B cells secreting antigen-specific IgA as well as B cells secreting IgA irrespective of antigen specificity (total IgA) can be measured.

Detection of antigen-specific B cells

Detection of antigen-specific IgA-secreting B cells, visualised as spots, can be made using either of two protocols (Protocol I or II). The choice of protocol will depend on several factors such as the available amount and molecular nature of the antigen. The use of biotinylated antigen (Protocol I) may be more sensitive and result in spots of higher quality. In addition, the antigen is used in smaller quantities and is not present during cultivation.

In **Protocol I**, antibodies from all secreting B cells are captured by anti-IgA antibodies coated on the ELISpot plate. Cells producing specific antibody are then visualised by the addition of biotinylated antigen followed by Streptavidin-enzyme and substrate. Please note that antigen must be biotinylated prior to starting the experiment for this protocol.

In **Protocol II**, representing the conventional way of performing B cell ELISpot, the antigen is coated on the ELISpot plate and antigen-specific antibodies secreted by the B cells bind to the immobilised antigen. Spots are detected after addition of biotinylated anti-IgA followed by Streptavidin-enzyme and substrate.

Detection of total IgA-producing B cells

Both protocols allow determination of the total number of IgA-secreting B cells, used as a positive control or for establishing the proportion of antigen-specific B cells.

Hints and comments

Please read through before starting the assay

These suggestions are based on the detection of antigen-specific immune responses using spleen cells. If using B-cell clones, mixtures of separated cell fractions etc., other protocols may have to be considered.

Plates

We recommend the use of PVDF-based membrane plates. Maximal antibody binding capacity of these plates is obtained by a brief treatment with ethanol. It is essential that the membrane is not allowed to dry after the treatment. If this occurs the treatment step (A2-A3) needs to be repeated before adding the coating antibody.

Plate washing

Washing of plates can be done using a multi-channel micropipette. In washing steps not requiring sterile conditions (C1-C5), a regular ELISA plate washer can also be used, provided that the washing head is adapted to the ELISpot plates. Avoid getting liquid on the underside of the membrane as this may cause leakage due to capillary drainage. If using MAIPSWU plates, always remove the plate from the plate tray before manually emptying the plate.

Cells

Both fresh and cryopreserved cells may be used with good results. However it is recommended that the latter are rested for at least one hour to allow removal of cell debris before addition to the plate. Triplicates or duplicates of 100,000-500,000 cells per well are often used to assess antigen-specific responses. For detection of total IgA, less cells (e.g. 50,000-100,000 cells per well) may be used to avoid confluent spot formation.

In vivo activated B cells

B cells that have been activated *in vivo*, for instance as a result of immunisation, may be analysed directly in the ELISpot wells without prior stimulation.

Memory B cells - *in vitro* activation

Memory B cells may require polyclonal stimulation before secreting detectable amounts of antibody. Pre-stimulation can be made with a mixture of R848 at 1 µg/ml and rmIL-2 at 10 ng/ml (both included) in separate tubes/plates for 48-72 hours. After pre-stimulation, wash the cells to ensure removal of any secreted antibodies. Resuspend the cells in medium and add the cell suspension to the ELISpot plate.

Reconstitution of rmIL-2: Add 1 ml PBS to obtain 0.5 µg/ml. Leave at room temperature for 15 minutes and then vortex. Use immediately or store in aliquots at -20°C.

Serum

The serum should be selected to support cell culture and give low background staining. We recommend the use of fetal calf serum. Alternatively serum-free medium evaluated for cell culture can be used.

Buffers

PBS for washing and dilution should be filtered (0.2 μm) for optimal results. Although possible to use, we do not recommend the inclusion of Tween or other detergents in the washing and incubation buffers.

Substrate development

Develop until distinct spots are visible in positive wells (usually 10-30 minutes). A general darkening of the membrane may occur but disappears after drying.

Biotinylation of antigen

If biotinylated antigen will be used for detection, (Protocol I), this should be labeled before starting the ELISpot experiment. A suitable reagent for biotinylation is e.g. Sulfo-NHS-LC-Biotin (Pierce product code 21335). Further information on biotin and biotin conjugation may be found at www.piercenet.com

-The protein antigen to be biotinylated with Sulfo-NHS-LC-Biotin should be in a suitable amine-free buffer (e.g. PBS, pH 7.4). If solution contains sodium azide, this should first be removed e.g. by dialysis. For small volumes (< 200 μl) a MINI-dialysis unit Slide-A-Lyser (10.000 MWCO, Pierce) can be used. The same dialysis unit can be used for removal of free biotin (see below). Conjugation is normally performed using a 20- to 50-fold molar excess of biotin to protein where the higher ratio is for protein at low concentration.

-The Sulfo-NHS-LC-Biotin should be dissolved in pure H_2O and immediately added to the antigen solution. After incubation for 30 minutes at room temperature, free biotin may be removed by dialysis against PBS or other buffer.

The concentration of the biotinylated antigen in the assay is typically 0.01-1 $\mu\text{g}/\text{ml}$ but depends on the size of the antigen and the degree of biotinylation and should be established separately.

Protocol I

(using biotinylated antigen)

A Preparation of ELISpot plate (sterile conditions)

1. Dilute the coating anti-IgA antibody to 10 µg/ml in sterile PBS, pH 7.4.
2. PVDF plates need to be treated (activated) with ethanol before coating. PVDF plates from Millipore, type MAIPSWU should be treated with 50 µl 70% ethanol/well for maximum 2 minutes. PVDF plates, type MSIP should be treated with 15 µl 35% ethanol/well for maximum 1 minute. After ethanol treatment, immediately proceed to plate washing.
3. Wash plate 5 times with sterile water, 200 µl/well.
4. Add 100 µl/well of the antibody solution and incubate overnight at 4-8°C.

B Incubation of cells in plate (sterile conditions)

See Hints and Comments section for details on cell stimulation.

1. Wash plate 5 times with sterile PBS, 200 µl/well, to remove excess antibody.
2. For blocking, add 200 µl/well of medium containing 10% of the same serum as used for the cell suspensions. Incubate for at least 30 minutes at room temperature.
3. For analysis of *in vivo* activated B cells, remove the medium and add the cell suspension to the ELISpot plate. For analysis of memory B cells, the cells can be pre-activated separately in tubes prior to addition to the plates. See further Hints and Comments/*In vivo* activated B cells and Memory B cells.
4. Put the plate in a 37°C humidified incubator with 5% CO₂ and incubate for 16-24 hours. Do not move the plate during this time and take measures to avoid evaporation (e.g. by wrapping the plate in aluminium foil).

C Detection of spots

1. To remove the cells, empty the plate and wash 5 times with PBS, 200 µl/well.
- 2a. **Antigen-specific IgA**
Dilute biotinylated antigen to a suitable concentration (0.01-1 µg/ml) in PBS-0.5% fetal calf serum (FCS). Add 100 µl/well and incubate for 2 hours at room temperature.
- 2b. **Total IgA**
Dilute the detection antibody (anti-IgA-biotin) to 1 µg/ml in PBS-0.5% FCS. Add 100 µl/well and incubate for 2 hours at room temperature.
3. Wash plate as above (step C1).
4. Dilute the Streptavidin-ALP (1:1000) in PBS-0.5% FCS and add 100 µl/well. Incubate for 1 hour at room temperature..
5. Wash plate as above (step C1).
6. Add 100 µl/well of substrate solution (e.g. BCIP/NBT) and develop until distinct spots emerge.
7. Stop colour development by washing extensively in tap water. If desirable, remove the plate from the plate tray or the underdrain and rinse the underside of the membrane.
8. Leave the plate to dry. Inspect and count spots in an ELISpot reader or in a dissection microscope.
9. Store plate in the dark at room temperature.

Protocol II

(using coated antigen)

A Preparation of ELISpot plate (sterile conditions)

- | | |
|--|---|
| 1a. Antigen-specific IgA
Dilute the antigen to suitable concentration (1-50 µg/ml in sterile PBS, pH 7.4). | 1b. Total IgA
Dilute the coating anti-IgA antibody to 10 µg/ml in sterile PBS, pH 7.4 |
|--|---|
2. PVDF plates need to be treated (activated) with ethanol before coating. PVDF plates from Millipore, type MAIPSWU should be treated with 50 µl 70% ethanol/well for maximum 2 minutes. PVDF plates, type MSIP should be treated with 15 µl 35% ethanol/well for maximum 1 minute. After ethanol treatment, immediately proceed to plate washing.
 3. Wash plate 5 times with sterile water, 200 µl/well.
 4. Add 100 µl/well of the antigen or antibody solution and incubate overnight at 4-8°C.

B Incubation of cells in plate (sterile conditions)

See Hints and Comments section for details on cell stimulation.

1. Wash plate 5 times with sterile PBS, 200 µl/well, to remove excess antibody
2. For blocking, add 200 µl/well of medium containing 10% of the same serum as used for the cell suspensions. Incubate for at least 30 minutes at room temperature.
3. For analysis of *in vivo* activated B cells, remove the medium and add the cell suspension to the ELISpot plate. For analysis of memory B cells, the cells can be pre-activated separately in tubes prior to addition to the plates. See further Hints and Comments/*In vivo* activated B cells and Memory B cells.
4. Put the plate in a 37°C humidified incubator with 5% CO₂ and incubate for 16-24 hours. Do not move the plate during this time and take measures to avoid evaporation (e.g. by wrapping the plate in aluminium foil).

C Detection of spots

1. To remove the cells, empty the plate and wash 5 times with PBS, 200 µl/well.
2. Dilute the detection antibody (anti-IgA-biotin) to 1 µg/ml in PBS-0.5% fetal calf serum (PBS-0.5% FCS). Add 100 µl/well and incubate for 2 hours at room temperature.
3. Wash plate as above (step C1).
4. Dilute the Streptavidin-ALP (1:1000) in PBS-0.5% FCS and add 100 µl/well. Incubate for 1 hour at room temperature.
5. Wash plate as above (step C1).
6. Add 100 µl/well of substrate solution (e.g. BCIP/NBT) and develop until distinct spots emerge.
7. Stop colour development by washing extensively in tap water. If desirable, remove the plate from the plate tray or the underdrain and rinse the underside of the membrane.
8. Leave the plate to dry. Inspect and count spots in an ELISpot reader or in a dissection microscope.
9. Store plate in the dark at room temperature.

MABTECH AB
Augustendalsvägen 19
Box 1233, SE-131 28 Nacka Strand
Sweden
Tel: +46 8 716 27 00
Fax: +46 8 716 27 01
E-mail: mabtech@mabtech.com

MABTECH Inc
M.E.B. 220, 3814 West Street
Cincinnati, OH 45227
USA
Toll free: +1 866 ELI-SPOT
Tel: +1 513 871 4500
Fax: +1 513 871 7353
E-mail: mabtech.usa@mabtech.com

MABTECH AUSTRALIA Pty Ltd
Australia
Tel: +61 3 9470 4704
Fax: +61 3 8678 3216
E-mail: mabtech.au@mabtech.com

MABTECH AB Büro Deutschland
Germany
Tel: +49 40 4135 7935
Fax: +49 40 4135 7945
E-mail: mabtech.de@mabtech.com

MABTECH AB Bureau de liaison
France
Tel: +33 (0)4 92 38 80 70
Fax: +33 (0)4 92 38 80 71
E-mail: mabtech.fr@mabtech.com

MABTECH

NOTE; for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages therefrom.



Intertek



Intertek

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the above standards.

For more information about Mabtech services and our products, visit

www.mabtech.com