

ELISpot Flex: Rabbit IgG (HRP)

Product Code: 3152-2H

CONTENTS:

Vial 1

Capture mAb: MT5131 (1.2 ml)

Concentration: 0.5 mg/ml

Vial 2

Detection mAb: MT536, biotin (50 µl)

Concentration: 0.5 mg/ml

Vial 3

Streptavidin-HRP (500 µl)


Antibodies are supplied in sterile filtered (0.2 µm) PBS with 0.02% sodium azide. Streptavidin-HRP is supplied in PBS with 0.002% Kathon CG. Vials have been overfilled to ensure recovery of the specified amount.

STORAGE:

Shipped at ambient temperature. On arrival antibodies and Streptavidin-HRP should be stored refrigerated at 4-8°C. The expiry date indicates how long unopened products, stored according to instructions, are recommended for use.


Schematic illustration of the different protocols

PROTOCOL I

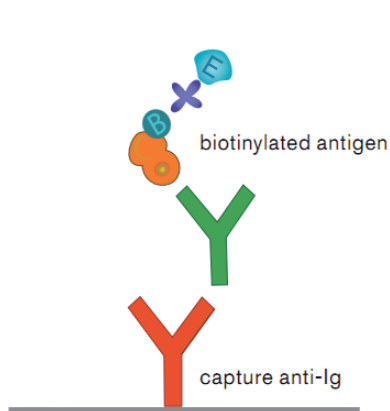

Antigen-specific Ig ELISpot

PROTOCOL II

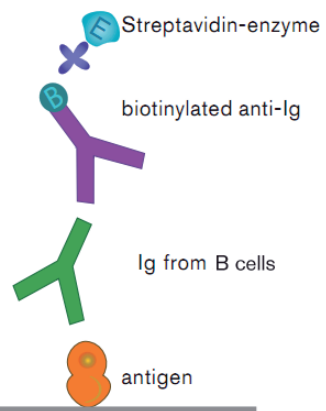

Total Ig ELISpot

Kit description

Please read through before starting the assay

This kit contains reagents for detection and enumeration of B cells secreting IgG antibodies. B cells secreting antigen-specific IgG as well as B cells secreting IgG irrespective of antigen specificity (total IgG) can be measured.

Detection of antigen-specific B cells

Detection of antigen-specific IgG-secreting B cells, visualised as spots, can be made using either of two protocols (Protocol I or II). The choice of protocol will depend on several factors such as the available amount and molecular nature of the antigen. The use of biotinylated antigen (Protocol I) may be more sensitive and result in spots of higher quality. In addition, the antigen is used in smaller quantities and is not present during cultivation.

Protocol I: Please note that antigen must be biotinylated prior to starting the experiment. IgG from all secreting B cells are captured by anti-IgG antibodies coated on the ELISpot plate. The secretion of specific IgG is then visualised by the addition of biotinylated antigen followed by Streptavidin-enzyme and substrate.

Protocol II: In the conventional way of performing B cell ELISpot, the antigen is coated on the ELISpot plate and antigen-specific IgG secreted by the B cells bind to the immobilised antigen. Spots are detected after addition of biotinylated anti-IgG followed by Streptavidin-enzyme and substrate.

Detection of total IgG-producing B cells

Both protocols allow determination of the total number of IgG-secreting B cells, used as a positive control or for establishing the proportion of antigen-specific B cells.

Hints and comments

Please read through before starting the assay

These suggestions are based on the detection of antigen-specific immune responses using PBMC and spleen cells. If using B-cell clones, mixtures of separated cell fractions etc., other protocols may have to be considered.

Plates

We recommend the use of PVDF-based membrane plates. Maximal antibody binding capacity of these plates is obtained by a brief treatment with ethanol. It is essential that the membrane is not allowed to dry after the treatment. If this occurs the treatment step (A2-A3) needs to be repeated before adding the coating antibody.

Plate washing

Washing of plates can be done using a multi-channel micropipette. In washing steps not requiring sterile conditions (C1-C5), a regular ELISA plate washer can also be used, provided that the washing head is adapted to the ELISpot plates. Avoid getting liquid on the underside of the membrane as this may cause leakage due to capillary drainage. If using MAIPSWU plates, always remove the plate from the plate tray before manually emptying the plate.

Cells

Both fresh and cryopreserved cells may be used with good results. However it is recommended that the latter are rested for at least one hour to allow removal of cell debris before addition to the plate. As a guideline, 100,000-500,000 cells per well may be used to assess antigen-specific responses. For detection of total IgG, less cells (e.g. 50,000-150,000 cells per well) may be used to avoid confluent spot formation. It is recommended to set up samples in triplicates or duplicates.

***In vivo* activated B cells**

The time-point for detection of cells secreting antigen-specific antibodies after immunization needs to be established by the user.

Memory B cells - *in vitro* activation

Serum

The serum should be selected to support cell culture and give low background staining. We recommend the use of fetal calf serum. Alternatively serum-free medium evaluated for cell culture can be used.

Buffers

PBS for washing and dilution should be filtered (0.2 µm) for optimal results. Although possible to use, we do not recommend the inclusion of Tween or other detergents in the washing and incubation buffers.

Substrate development

Develop until distinct spots are visible in positive wells (usually 5-30 minutes). A general darkening of the membrane may occur but disappears after drying. Preferably use deionized water to stop the plates since some ions may cause fading of TMB spots.

Biotinylation of antigen

If biotinylated antigen will be used for detection, (Protocol I), this should be labeled before starting the ELISpot experiment. A suitable reagent for biotinylation is e.g. Sulfo-NHS-LC-Biotin (Pierce product code 21335). Further information on biotin and biotin conjugation may be found at www.piercenet.com

-The protein antigen to be biotinylated with Sulfo-NHS-LC-Biotin should be in a suitable amine-free buffer (e.g. PBS, pH 7.4). If solution contains sodium azide, this should first be removed e.g. by dialysis. For small volumes (< 200 µl) a MINI-dialysis unit Slide-A-Lyser (10.000 MWCO, Pierce) can be used. The same dialysis unit can be used for removal of free biotin (see below). Conjugation is normally performed using a 20- to 50-fold molar excess of biotin to protein where the higher ratio is for protein at low concentration.

-The Sulfo-NHS-LC-Biotin should be dissolved in pure H₂O and immediately added to the antigen solution. After incubation for 30 minutes at room temperature, free biotin may be removed by dialysis against PBS or other buffer.

The concentration of the biotinylated antigen in the assay is typically 0.01-1 µg/ml but depends on the size of the antigen and the degree of biotinylation and should be established separately.

Protocol I

(using biotinylated antigen)

A Preparation of ELISpot plate (sterile conditions)

1. Dilute the coating mAb MT5131 to 15 µg/ml in sterile PBS, pH 7.4.
2. PVDF plates need to be treated (activated) with ethanol before coating. PVDF plates from Millipore, type MAIPSWU should be treated with 50 µl 70% ethanol/well for maximum 2 minutes. PVDF plates, type MSIP should be treated with 15 µl 35% ethanol/well for maximum 1 minute. After ethanol treatment, immediately proceed to plate washing.
3. Wash plate 5 times with sterile water, 200 µl/well.
4. Add 100 µl/well of the antibody solution and incubate overnight at 4-8°C.

B Incubation of cells in plate (sterile conditions)

See Hints and Comments section for details on cell stimulation.

1. Wash plate 5 times with sterile PBS, 200 µl/well, to remove excess antibody.
2. For blocking, add 200 µl/well of medium containing 10% of the same serum as used for the cell suspensions. Incubate for at least 30 minutes at room temperature.
3. Remove the medium and add the cell suspension to the ELISpot plate (100 µl/well).
4. Put the plate in a 37°C humidified incubator with 5% CO₂ and incubate for 16-24 hours. Do not move the plate during this time and take measures to avoid evaporation (e.g. by wrapping the plate in aluminium foil).

C Detection of spots

1. To remove the cells, empty the plate and wash 5 times with PBS, 200 µl/well.
- 2a. **Antigen-specific IgG**
Dilute biotinylated antigen to a suitable concentration (0.01-1 µg/ml) in PBS-0.5% fetal calf serum (FCS). Add 100 µl/well and incubate for 2 hours at room temperature.
- 2b. **Total IgG**
Dilute the detection mAb MT536-biotin to 0.5 µg/ml in PBS-0.5% FCS. Add 100 µl/well and incubate for 2 hours at room temperature.
3. Wash plate as above (step C1).
4. Dilute the Streptavidin-HRP* in PBS-0.5% FCS and add 100 µl/well. Incubate for 1 hour at room temperature.
*The dilution of Streptavidin-HRP should be adapted to the substrate, e.g., 1:1000 for TMB and 1:100 for AEC. *Please note that sodium azide used in buffers will inhibit HRP activity.*
5. Wash plate as above (step C1).
6. Add 100 µl/well of substrate solution (e.g. TMB) and develop until distinct spots emerge.
7. Stop colour development by washing extensively in deionized water. If desirable, remove the plate from the plate tray or the underdrain and rinse the underside of the membrane.
8. Leave the plate to dry. Inspect and count spots in an ELISpot reader or in a dissection microscope.
9. Store plate in the dark at room temperature.

Protocol II

(using coated antigen)

A Preparation of ELISpot plate (sterile conditions)

- | | |
|--|--|
| 1a. Antigen-specific IgG
Dilute the antigen to suitable concentration (1-50 µg/ml in sterile PBS, pH 7.4). | 1b. Total IgG
Dilute the coating mAb MT5131 to 15 µg/ml in sterile PBS, pH 7.4 |
|--|--|
2. PVDF plates need to be treated (activated) with ethanol before coating. PVDF plates from Millipore, type MAIPSWU should be treated with 50 µl 70% ethanol/well for maximum 2 minutes. PVDF plates, type MSIP should be treated with 15 µl 35% ethanol/well for maximum 1 minute. After ethanol treatment, immediately proceed to plate washing.
 3. Wash plate 5 times with sterile water, 200 µl/well.
 4. Add 100 µl/well of the antigen or antibody solution and incubate overnight at 4-8°C.

B Incubation of cells in plate (sterile conditions)

See Hints and Comments section for details on cell stimulation.

1. Wash plate 5 times with sterile PBS, 200 µl/well, to remove excess antibody
2. For blocking, add 200 µl/well of medium containing 10% of the same serum as used for the cell suspensions. Incubate for at least 30 minutes at room temperature.
3. Remove the medium and add the cell suspension to the ELISpot plate (100 µl/well).
4. Put the plate in a 37°C humidified incubator with 5% CO₂ and incubate for 16-24 hours. Do not move the plate during this time and take measures to avoid evaporation (e.g. by wrapping the plate in aluminium foil).

C Detection of spots

1. To remove the cells, empty the plate and wash 5 times with PBS, 200 µl/well.
2. Dilute the detection mAb MT536 to 0.5 µg/ml in PBS-0.5% fetal calf serum (PBS-0.5% FCS). Add 100 µl/well and incubate for 2 hours at room temperature.
3. Wash plate as above (step C1).
4. Dilute the Streptavidin-HRP* in PBS-0.5% FCS and add 100 µl/well. Incubate for 1 hour at room temperature.
*The dilution of Streptavidin-HRP should be adapted to the substrate, e.g., 1:1000 for TMB and 1:100 for AEC. *Please note that sodium azide used in buffers will inhibit HRP activity.*
5. Wash plate as above (step C1).
6. Add 100 µl/well of substrate solution (e.g. TMB) and develop until distinct spots emerge.
7. Stop colour development by washing extensively in deionized water. If desirable, remove the plate from the plate tray or the underdrain and rinse the underside of the membrane.
8. Leave the plate to dry. Inspect and count spots in an ELISpot reader or in a dissection microscope.
9. Store plate in the dark at room temperature.

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the standards ISO 9001:2015 & ISO 13485:2016.



The products are for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages there from.

Mabtech AB (Head Office)

Sweden

Tel: +46 8 716 27 00

mabtech@mabtech.com

Mabtech, Inc.

USA

Tel: +1 513 871-4500

mabtech.usa@mabtech.com

www.mabtech.com