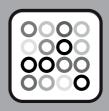


ELISA Pro: Human Latent TGF-β1

3550-1HP-1 | 3550-1HP-2 | 3550-1HP-10

Datasheet & Protocol






Table of Contents

| Introduction | . 4 |
|----------------------|-----|
| Shipping and storage | . 4 |
| Contents | . 5 |
| Preparation | . 6 |
| Protocol | . 8 |
| Performance | 10 |

Introduction

Mabtech's carefully validated ELISA Pro kits provide all the necessary reagents to conveniently quantify analytes in serum, plasma, and cell culture supernatants in a robust, sensitive, and specific manner. For Research Use Only (RUO). Not for use in diagnostic procedures.

ELISA assay principle

ELISA Pro kits are supplied with ELISA strip plates precoated with monoclonal antibody (mAb). Analyte in the sample is captured by the coated mAb and detected by the biotinylated detection mAb followed by Streptavidin-HRP (SA-HRP). Addition of TMB substrate will result in a colored substrate product. The reaction is stopped with sulfuric acid and the optical density can be quantified using an ELISA plate reader. The concentration of analyte is determined by comparison to a serial dilution of the ELISA standard analyzed in parallel.

Analysis of serum and plasma samples

The ELISA Pro kits include ELISA diluent, a buffer that prevents false-positive signals. The buffer blocks heterophilic antibodies from cross-linking the assay antibodies. Heterophilic antibodies are commonly found

les

Indicate the latest three three

in human serum/plasma and can also be present in other species. The buffer has been validated using serum/plasma samples from healthy human blood donors.

Shipping and storage

The kit is shipped at ambient temperature. All reagents should be stored at 4-8 °C upon receipt, except the standard, which should be stored at -20 °C. The expiry date indicates how long unopened products, stored according to instructions, are recommended for use. Do not combine components from different kit batches or components from other suppliers.

Contents

| Component | 1-plate kit | 2-plate kit | 10-plate kit | |
|--|--------------|--------------|---------------|--|
| Pre-coated ELISA strip plate: Anti-TGF-β1 mAb MT593 | 1 x 96 wells | 2 x 96 wells | 10 x 96 wells | |
| Recombinant human LAP ELISA standard | 1 vial | 1 vial | 1 vial | |
| Detection mAb MT517, biotin (0.5 mg/ml) | 25 μΙ | 50 μΙ | 250 μΙ | |
| Streptavidin-HRP | 15 μΙ | 25 μΙ | 125 μΙ | |
| Standard reconstitution buffer A8 | 1 ml | 1 ml | 1 ml | |
| Wash buffer concentrate | 120 ml | 120 ml | 5 x 120 ml | |
| ELISA diluent | 120 ml | 120 ml | 3 x 120 ml | |
| Streptavidin-HRP diluent | 15 ml | 25 ml | 120 ml | |
| TMB substrate | 15 ml | 25 ml | 120 ml | |
| Stop solution | 15 ml | 25 ml | 120 ml | |
| Adhesive plate covers | 3 | 6 | 30 | |

To ensure total recovery of the stated quantity, bottles and vials have been overfilled.

Materials required but not supplied

- Microplate reader capable of reading at 450 nm
- ELISA plate washer; automated or manual (e.g., multipipette or squirt bottle)
- · Precision pipettes, tips, and graduated cylinders
- Tubes for standard and sample dilutions
- · Distilled or deionized water

Safety information

The Stop solution, 0.18 M ${\rm H_2SO_4}$ (< 1%), is irritating to eyes and skin and should be handled with care. The standard should also be handled carefully as the effects of exposure are unknown. Buffers and reagents in solution contain the preservative Kathon CG (0.002%), a potential allergen that may cause sensitization through skin contact. Human and animal samples should be treated as potentially hazardous biologic material. All material should be disposed of in accordance with local regulations. For further information please consult the Safety Data Sheet on our website.

Preparation

- Allow the plates and assay reagents to reach room temperature before starting the assay (except for the TMB substrate which should preferably be used cold).
- Plan the plate layout to include a standard curve, samples, and an assay background control, all in duplicate. The volume per well should not exceed 100 µl. Include a plate blank (wells with only Substrate and Stop solution) to be used for subtraction before analysis.

Wash buffer

Add 50 ml Wash buffer concentrate to 950 ml distilled or deionized water (sufficient for all washing steps of 1 plate). If crystals have formed in the 20x concentrate, bring to room temperature and mix gently to dissolve.

Samples


For accurate assay performance, all samples should be diluted at least 2-fold in ELISA diluent. Samples can be diluted in tubes or plates; diluent should be added before the samples and visible precipitates should be removed. The use of strongly hemolyzed and hyperlipemic samples may give inaccurate determination of the concentration. Samples containing high levels of analyte exceeding the standard range of the assay will require further dilution. For quantification of latent TGF- β 1 in blood, the use of plasma is recommended since serum contains high levels of latent TGF- β 1 released from platelets during sample preparation. Plasma can be obtained using EDTA, citrate or heparin as anticoagulants. To minimize the platelet content in the sample, and additional centrifugation of plasma at 10,000 x g for 10 minutes is recommended.

ELISA standard

Recombinant human LAP (homodimer) is used as standard. Since it differs in molecular weight from latent TGF- β 1. determination of latent TGF- β 1 using the standard curve is based on a molar comparison. Reconstitute the standard to a stock solution of 20 nM by adding 1 ml of Standard reconstitution buffer. Allow the standard to dissolve for 5 minutes and mix thoroughly. The standard should be kept in aliquots at -20 °C. Avoid repeated freeze-thaw cycles.

Preparation of standard curve

Dilute the standard stock solution to create a standard curve as shown. The indicated volumes are sufficient for duplicates. The last vial is used as an assay background control, i.e., the standard should be omitted. Prepare the standard curve within 30 minutes of use.

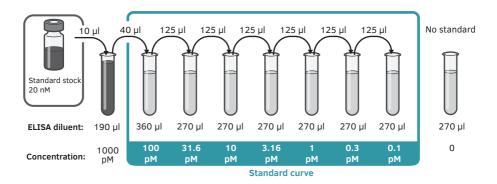

Detection antibody

Dilute the detection mAb in ELISA diluent to 1 μ g/ml within 15 minutes of use. For each plate, add 24 μ l detection mAb to 12 ml ELISA diluent.

Streptavidin-HRP


Dilute the Streptavidin-HRP 1:1000 in Streptavidin-HRP diluent within 15 minutes of use. For each plate, add 12 μ l Streptavidin-HRP to 12 ml Streptavidin-HRP diluent.

Protocol

Prepare the reagents, standard curve, and samples as described in the Preparation section. Assemble the required number of strips in the plate frame and label the top of each strip. Store the remaining strips in the foil bag containing the desiccant at 4-8 °C.


- 1. Wash the plate 5 times with wash buffer, 300 μ l per well. After the final wash, invert and tap the plate firmly against absorbent paper. Immediately proceed to the next step.
- **2.** Add 100 µl per well of samples (diluted at least 2-fold), standard, and assay background control. Mix by tapping the plate. Cover the plate with an adhesive plate cover and incubate for 2 hours at room temperature.
- **3.** Wash as in step 1.
- **4.** Add 100 µl per well of detection mAb. Cover the plate and incubate for 1 hour at room temperature.
- **5.** Wash as in step 1.
- **6.** Add 100 µl per well of Streptavidin-HRP. Cover the plate and incubate for 1 hour at room temperature.
- **7.** Wash as in step 1.
- **8.** Add 100 µl of TMB substrate to each well. Incubate at room temperature, protected from direct light for 15 minutes.
- **9.** Add 100 μ l of Stop solution to each well to stop the color development.
- **10.** Measure absorbance at 450 nm within 15 minutes. Preferably use a reader capable of subtracting a reference wavelength between 570 and 650 nm.

We recommend the use of an ELISA software utilizing a 4- or 5-parameter curve fit. Subtract the mean absorbance value of the blank from the samples, standard and assay background control prior to creating the standard curve and analyzing the results.

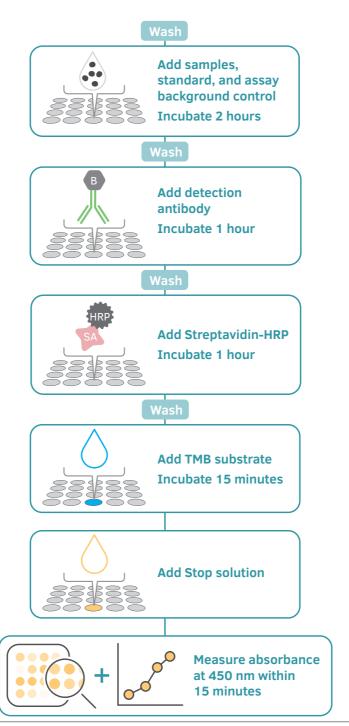

Performance



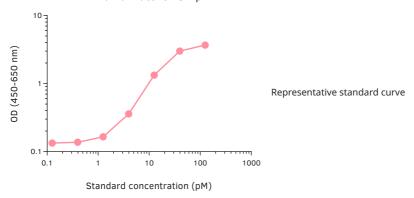

Standard range 0.4-40 pM

Sensitivity 0.13 pM

The limit of detection is the lowest concentration that is possible to detect but not necessarily quantify with precision and accuracy.

Calibration

Recombinant human LAP (homodimer) is used as standard. Since it differs in molecular weight from latent TGF- β 1, determination of latent TGF- β 1 using the standard curve is based on a molar comparison. A concentration of 1 pM LAP corresponds to 1 pM latent TGF- β 1. For conversion from pM to pg/ml: 1 pM LAP = 54 pg/ml and 1 pM latent TGF- β 1 = 80 pg/ml. Note that cytokine values obtained should be multiplied with the dilution factor used for each sample

Precision

| | Intra-assay | | | Inter-assay | | |
|-----------|-------------|------|------|-------------|------|------|
| Sample | 1 | 2 | 3 | 1 | 2 | 3 |
| n | 8 | 8 | 8 | 3 | 3 | 3 |
| Mean (pM) | 24.6 | 11.3 | 6.7 | 23.9 | 12.8 | 7.6 |
| SD | 0.7 | 0.3 | 0.07 | 2.0 | 1.3 | 0.8 |
| CV% | 3.0 | 2.2 | 1.1 | 8.4 | 10.4 | 10.9 |

Intra-assay and inter-assay precision were determined at 3 different concentrations of analyte (8 replicates per concentration in 3 assays).

Linearity

Dilution of a human plasma gives a mid-curve recovery of 89-96% in repeated experiments (mean 94%).

Specificity

TGF- β 1 is secreted in a latent, biologically inactive form, non-covalently linked with a Latency Associated Protein (LAP). The mAbs MT593 and MT517-biotin bind to the LAP entity of human latent TGF- β 1. This enables direct quantification of latent TGF- β 1 in ELISA with no need of sample pre-treatment. The ELISA does not detect human latent TGF- β 2 or - β 3 or bovine latent TGF- β 1.

The mAbs in this kit cross-react with latent TGF-β1 from non-human primates. Please visit www.mabtech.com for reactivity on NHP species.



Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the standards ISO 9001:2015 & ISO 13485:2016.





The products are for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages there from.

Mabtech AB (Head Office) Sweden

Sweden Tel: +46 8 716 27 00 mabtech@mabtech.com Mabtech, Inc.

Tel: +1 513 871-4500 mabtech.usa@mabtech.com