

FluoroSpot Plus: Human IL-6/IL-1 β /GM-CSF/TNF- α

FSP-13121509-2
FSP-13121509-10

www.mabtech.com

MABTECH

The products are for research use only.


MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages therefrom.

INTRODUCTION

FluoroSpot - General description

This four analyte FluoroSpot kit enables convenient enumeration of cells secreting either one, two, three or four analytes.

FluoroSpot combines the sensitivity of ELISpot with the capacity to analyze secretion of several analytes simultaneously. This highly sensitive cellular assay is robust, easy to perform, and suitable both for testing single samples and for large-scale screening.

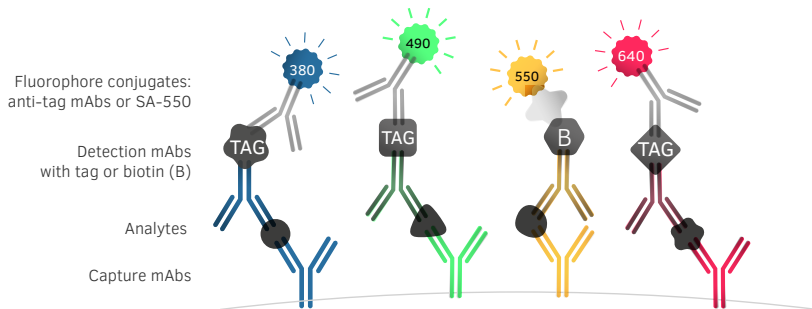

Our FluoroSpot Plus kits include pre-coated plates ready to use. Monoclonal antibodies (mAbs) are coated onto a PVDF membrane in a 96-well plate. If several analytes are analyzed, mAbs with different specificities are coated. Detection of cells secreting single or multiple analytes is made possible using biotinylated and/or tag-labeled detection mAbs. The detection step is visualized and amplified by specific fluorophore-conjugated reagents and the resulting spots are analyzed in an automated reader.

CONTENTS AND STORAGE

FluoroSpot Plus: Human IL-6/IL-1 β /GM-CSF/TNF- α

The kit is shipped at ambient temperature. All reagents should be stored at 4-8°C upon arrival. Plates should be stored at room temperature. The expiry date indicates how long unopened products, stored according to instructions, are recommended to remain in use.

	Component	Concentration	2-plate kit	10-plate kit
Detection mAbs	Anti-IL-6 mAb (39C3), DIG		120 μ l	600 μ l
	Anti-IL-1 β mAb (7P10), BAM		120 μ l	600 μ l
	Anti-GM-CSF mAb (23B6), biotin	1 mg/ml	50 μ l	250 μ l
	Anti-TNF- α mAb (MT20D9), WASP		120 μ l	600 μ l
Fluorophore conjugates	Anti-DIG mAb, 380		120 μ l	600 μ l
	Anti-BAM mAb, 490		120 μ l	600 μ l
	SA-550		120 μ l	600 μ l
	Anti-WASP mAb, 640		120 μ l	600 μ l
Stimulus	LPS	0.05 mg/ml	120 μ l	120 μ l
Buffer/Solution	FluoroSpot enhancer		25 ml	60 ml
Pre-coated plates	Plates pre-coated with: mAbs 13A5, MT175, 21C11 and MT25C5		2	10

The biotinylated antibody contains 0.02% sodium azide. The DIG-, BAM-, and WASP-conjugated antibodies, fluorophore-conjugated reagents, and the FluoroSpot enhancer contain 0.002 % Kathon CG. Vials have been overfilled to ensure recovery of stated volume.

ASSAY PROTOCOL

Please consult the table below for recommendations on dilutions and cell incubation time. The reagents should not be diluted in advance. For additional information, read the **Tips For A Successful Assay** section.

	Component	Dilution	Step
Detection mAbs	Anti-IL-6 mAb (39C3), DIG	1:200	C2
	Anti-IL-1 β mAb (7P10), BAM	1:200	C2
	Anti-GM-CSF mAb (23B6), biotin	1:500	C2
	Anti-TNF- α mAb (MT20D9), WASP	1:200	C2
Fluorophore conjugates	Anti-DIG mAb, 380	1:200	C4
	Anti-BAM mAb, 490	1:200	C4
	SA-550	1:200	C4
	Anti-WASP mAb, 640	1:200	C4
Stimulus	LPS	1:1000	B3
Cell incubation time overnight			B2

A. Plate preparation (sterile conditions)

1. Remove the plate from the sealed package and wash the plate three times with sterile PBS (200 μ l/well).
2. Block/condition the plate by adding cell incubation medium containing 10% fetal calf serum (200 μ l/well). Incubate for at least 30 minutes at room temperature.

B. Cell incubation (sterile conditions)

1. Remove the medium and add the stimuli followed by the cell suspension. Cells and stimuli can also be mixed before addition to the plate.
2. Incubate the cells in the plate at 37 $^{\circ}$ C in a humidified incubator with 5% CO₂. Approximate incubation time is indicated above and can be further optimized by the user. Do not move the plate during incubation and avoid evaporation, for example by wrapping the plate in aluminum foil.

C. Detection

1. Remove the cells by emptying the plate and then wash the plate five times with PBS (200µl/well).
2. Dilute the detection antibodies in PBS containing 0.1% BSA (PBS-0.1% BSA) according to the table above. If more than one analyte is analyzed, dilute the detection antibodies in the same tube. Add 100 µl/well and incubate for 2 hours at room temperature.
3. Wash the plate five times with PBS (200 µl/well).
4. Dilute the fluorophore conjugates in PBS-0.1% BSA to the concentrations as shown in the table above. If more than one analyte is analyzed, dilute the conjugates in the same tube. Add 100 µl/well and incubate for 1 hour at room temperature. Protect the plate from light throughout the assay.
5. Wash the plate five times with PBS (200 µl/well).
6. Empty the plate, add FluoroSpot enhancer (50 µl/well), and leave for 5-15 minutes at room temperature.
7. Empty by flicking the plate to remove the FluoroSpot enhancer. Do not wash. It is important that the plate is not tapped against paper towels since dust particles may interfere with the assay.
8. Remove the underdrain (the soft plastic drain under the plate), for example using pliers. Dry the plate protected from light. The plate should be completely dry before analysis. Store plate in the dark at room temperature.
9. Spot analysis is performed with an automated FluoroSpot reader equipped with filters for the fluorophores used. Filters should have high specificity to avoid bleed-through artifacts.

Analyte	Required filter
IL-6	Excitation 380 nm/emission 430 nm (LED380)
IL-1 β	Excitation 490 nm/emission 510 nm (LED490)
GM-CSF	Excitation 550 nm/emission 570 nm (LED550)
TNF- α	Excitation 640 nm/emission 660 nm (LED640)

TIPS FOR A SUCCESSFUL ASSAY

Buffers

PBS for washing and dilution should be filtered (0.2 µm) to remove any particles. We do not recommend the inclusion of Tween or other detergents in the washing and incubation buffers.

Washing

Plates can be washed using a multi-channel micropipette. A regular plate washer can be used in washing steps when sterile conditions are not required, provided that the washing head is adapted to ELISpot/FluoroSpot plates.

Serum

Serum should be selected to support cell culture and give low background staining. We recommend the use of fetal calf serum. Alternatively, serum-free medium evaluated for cell culture can be used.

Cells

Both freshly prepared and cryopreserved cells may be used in FluoroSpot. We recommend that cryopreserved cells are rested for at least 1 hour after thawing to allow removal of cell debris before addition to the plate. Triplicates or duplicates of 200,000-500,000 cells/well are usually recommended to assess antigen-specific responses. For LPS-stimulated human PBMC, around 5000 cells/well is recommended. For purified LPS-stimulated monocytes, fewer cells will be required to avoid confluent spot formation.


Assay controls

The number of cells responding to stimulation is often compared to the number of cells spontaneously producing the cytokine, which is determined by incubating the same number of cells in the absence of stimuli. A polyclonal activator such as LPS (1-10 µg/ml) is often used as a control for cell viability and functionality of the test system.

Analysis

An automated FluoroSpot reader with separate filters is required for analysis. The reader should be equipped with filters for the fluorophores used. Filters should be selective for the specific wavelengths to avoid bleed-through artifacts. Spots from cells secreting multiple analytes are identified by coordinates, i.e. based on spot position. In order to obtain accurate measurements of cells secreting multiple analytes, the automated reader must be able to create an overlay analysis from the filters. To avoid fading, the plates should be kept dark and should preferably be analyzed within 1 week of development.

PROTOCOL SUMMARY

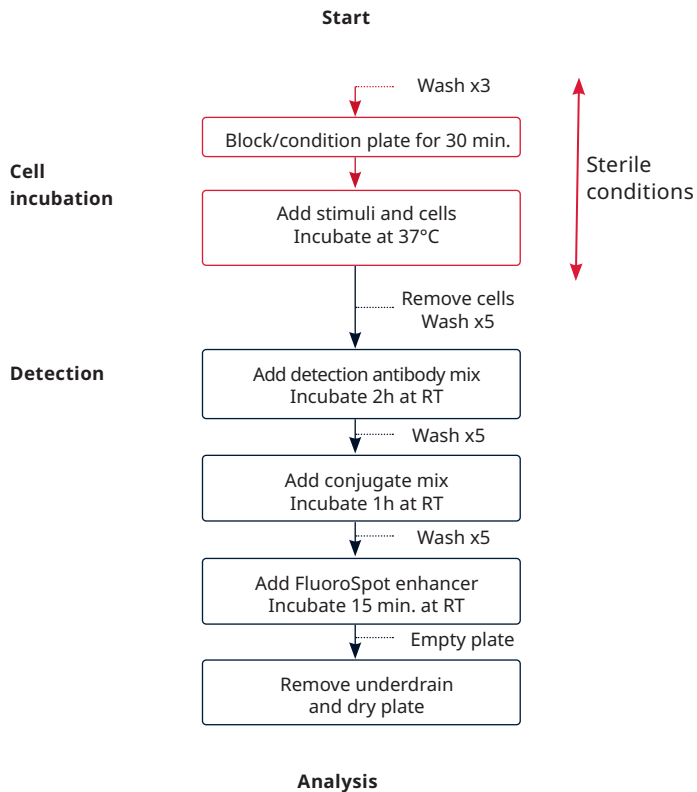

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the standards ISO 9001:2015 & ISO 13485:2016.



The products are for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages there from.

Mabtech AB (Head Office)
Sweden
Tel: +46 8 716 27 00
mabtech@mabtech.com

Mabtech, Inc.
USA
Tel: +1 513 871-4500
mabtech.usa@mabtech.com