

# Ferret IFN- $\gamma$ ELISA development kit

Product Code: 31 12-1H-20

---

CONTENTS, development kit for 20 plates:

**Vial 1 (green top)**

Monoclonal antibody MTF14 (1 ml)

Concentration: 0.5 mg/ml

**Vial 2 (yellow top)**

Biotinylated monoclonal antibody MTF19 (500  $\mu$ l)

Concentration: 0.5 mg/ml

**Vial 3 (white top)**

Streptavidin-Horseradish Peroxidase (250  $\mu$ l)

**Vial 4**

Recombinant ferret IFN- $\gamma$  standard

To ensure total recovery of stated quantity, vials have been overfilled.

**STORAGE:**

Shipped at ambient temperature. On arrival box 1 should be stored refrigerated at 4-8°C and box 2 should be stored frozen at -20°C.


# General

**Intended use:** For quantitative determination of native and recombinant ferret IFN- $\gamma$  in solution, e.g. cell culture supernatant and serum/plasma samples.

**Reagents:** Antibodies are supplied in sterile-filtered (0.2  $\mu\text{m}$ ) PBS with sodium azide (0.02%). Streptavidin-HRP is supplied in PBS with 1% BSA and 0.002% Kathon CG.

**Standard range:** 12-1200 pg/ml.

**Standard calibration:** No international standard exists for calibration

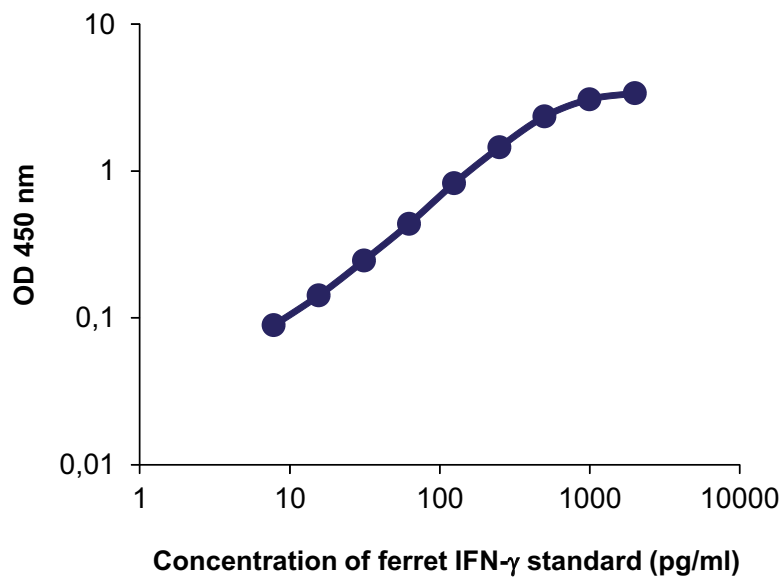

# Guidelines for Ferret IFN- $\gamma$ ELISA

- Day 1**
1. Coat a high protein binding ELISA plate with mAb MTF14, diluted to 2  $\mu\text{g}/\text{ml}$  in PBS, pH 7.4, by adding 100  $\mu\text{l}/\text{well}$ . Incubate overnight at 4-8°C.
- Day 2**
2. Wash twice with PBS (200  $\mu\text{l}/\text{well}$ ).
  3. Block plate by adding 200  $\mu\text{l}/\text{well}$  of PBS with 0.05% Tween 20 containing 0.1% BSA (incubation buffer). Incubate for 1 hour at room temperature.
  4. Wash five times with PBS containing 0.05% Tween.
  5. Prepare ferret IFN- $\gamma$  standard by reconstituting contents of vial 4 in 1 ml PBS with 1% BSA and leave at room temperature for 15 minutes, then vortex the tube and spin down. This gives a concentration of 0.5  $\mu\text{g}/\text{ml}$ . Use immediately or store in aliquots at -20°C for future use. We recommend the aliquots not to be refrozen after initial use. For the test, prepare dilutions of the stock using the standard range as a guideline.
  6. Add 100  $\mu\text{l}/\text{well}$  of samples or standards diluted in incubation buffer and incubate for 2 hours at room temperature.
  7. Wash as in step 4.
  8. Add 100  $\mu\text{l}/\text{well}$  of mAb MTF19-biotin at 1  $\mu\text{g}/\text{ml}$  in incubation buffer. Incubate for 1 hour at room temperature.
  9. Wash as in step 4.
  10. Add 100  $\mu\text{l}/\text{well}$  of Streptavidin-HRP diluted 1:1000 in incubation buffer. Incubate for 1 hour at room temperature. **Please note that sodium azide used in buffers will inhibit HRP activity.**
  11. Wash as in step 4.
  12. Add 100  $\mu\text{l}/\text{well}$  of appropriate substrate solution e.g. TMB, available from Mabtech, product code 3652-F10.
  13. Measure the optical density in an ELISA reader after suitable developing time. If required stop the reaction first.

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the standards ISO 9001:2015 & ISO 13485:2016.



**The products are for research use only.**

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages there from.

Mabtech AB (Head Office)  
Sweden  
Tel: +46 8 716 27 00  
mabtech@mabtech.com

Mabtech, Inc.  
USA  
Tel: +1 513 871-4500  
mabtech.usa@mabtech.com