

**Human Interferon- γ
Monoclonal antibody**
FOR NEUTRALIZATION

General

Antibody:	MT111W
Product code:	3420-1N-500
Quantity:	500 µg
Immunogen:	Recombinant human IFN- γ
Isotype:	Rat IgG2a
Specificity:	Native and recombinant human IFN- γ . Cross-reacts with IFN- γ from non-human primates (NHP). Please visit www.mabtech.com for detailed information on NHP reactivity.
Concentration:	1 mg/ml
Buffer:	Supplied in sterile filtered (0.2 µm) PBS. Contains no preservatives. The endotoxin content in the product is below 10 EU/mg.
Purification:	Purified from <i>in vitro</i> cultures by protein G affinity chromatography.
Storage:	Store product at 4-8°C or frozen at -20°C or below. Avoid repeated freezing/thawing.
Applications:	Neutralization of human IFN- γ bioactivity. It is recommended to establish the optimal antibody concentration for the assay system used. The antibody is also suitable for Western blot.


Reference: Epitope Mapping of Neutralizing Monoclonal Antibodies to Human Interferon- γ Using Human-Bovine Interferon- γ Chimeras. Zuber B. et.al. J. Interferon Cytokine Res. 2016 Sep;36(9):542-51.

Neutralization

Neutralization assay: Neutralization of IFN- γ induced activation of HEK BLUE IFN- γ sensor cells. To determine the neutralizing capacity and ND₅₀, HEK BLUE IFN- γ sensor cells were incubated with pre-mixtures of recombinant human IFN- γ (100 pg/ml) with or without serial dilutions of mAb MT111W. IFN- γ -induced enzyme secretion was measured after 20 h.

Neutralization dose 50 (ND₅₀): MT111W neutralized the effect of IFN- γ at a ND₅₀ value of 0.004-0.013 μ g/ml when 100 pg/ml IFN- γ was present.

ND₅₀ is defined as the concentration of antibody that results in 50% inhibition of the biological response under the experimental conditions specified.

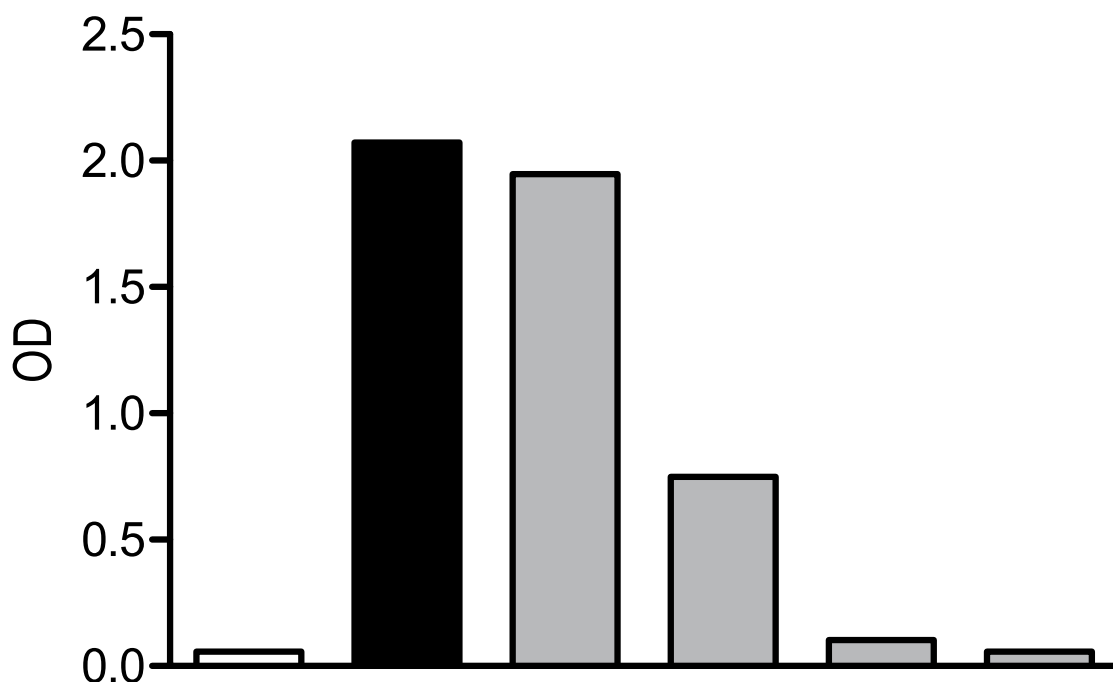

IFN- γ (100 pg/ml)

-

+

+

+

+

+

MT111W (μ g/ml)

-

-

0.002

0.02

0.2

2

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the standards ISO 9001:2015 & ISO 13485:2016.



The products are for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages there from.

Mabtech AB (Head Office)
Sweden
Tel: +46 8 716 27 00
mabtech@mabtech.com

Mabtech, Inc.
USA
Tel: +1 513 871-4500
mabtech.usa@mabtech.com