Mouse IL-10 ELISA development kit

Product Code: 3432-1A-20

CONTENTS, development kit for 20 plates:

Vial 1 (green top)

Monoclonal antibody MT60 (1000 µl)

Concentration: 0.5 mg/ml

Vial 2 (yellow top)

Biotinylated monoclonal antibody 51E6 (250 µl)

Concentration: 0.5 mg/ml

Vial 3 (white top)

Streptavidin-Alkaline Phosphatase (250 µl)

Vial 4

Recombinant mouse IL-10 standard

To ensure total recovery of stated quantity, vials have been overfilled.

STORAGE:

Shipped at ambient temperature. On arrival box 1 should be stored refrigerated at 4-8°C and box 2 should be stored frozen at -20°C.


General

Intended use: For quantitative determination of native and recombinant mouse IL-10 in solution, e.g. cell culture supernatant.

Reagents: Antibodies are supplied in sterile-filtered (0.2 μ m) PBS with sodium azide (0.02%). Streptavidin-ALP is supplied in 0.1 M Tris buffer with 0.002% Kathon CG.

Standard range: 4-400 pg/ml

Standard calibration: No international standard exists for calibration.

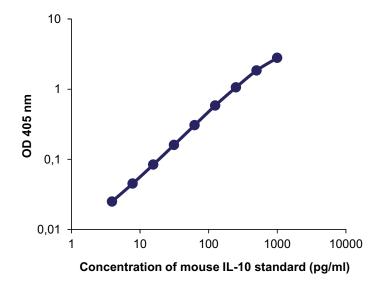

Guidelines for Mouse IL-10 ELISA

- **Day 1** 1. Coat a high protein binding ELISA plate with mAb MT60, diluted to 2 μg/ml in PBS, pH 7.4, by adding 100 μl/well. Incubate overnight at 4-8°C.
- **Day 2** 2. Wash twice with PBS (200 μl/well).
 - 3. Block plate by adding 200 μ l/well of PBS with 0.05% Tween 20 containing 0.1% BSA (incubation buffer). Incubate for 1 hour at room temperature.
 - 4. Wash five times with PBS containing 0.05% Tween.
 - 5. Prepare mouse IL-10 standard by reconstituting contents of vial 4 in 1 ml PBS with 1% BSA to give a concentration of 0.5 μg/ml. Leave at room temperature for 15 minutes and then vortex the tube and spin down. Use immediately or store in aliquots at -20°C for future use. We recommend the aliquots not to be refrozen after initial use. For the test, prepare dilutions of the stock using the standard range as a guideline.
 - 6. Add 100 μl/well of samples or standards diluted in incubation buffer and incubate for 2 hours at room temperature.
 - 7. Wash as in step 4.
 - 8. Add 100 μ l/well of mAb 51E6-biotin at 0.5 μ g/ml in incubation buffer. Incubate for 1 hour at room temperature.
 - 9. Wash as in step 4.
 - 10. Add 100 μl/well of Streptavidin-ALP diluted 1:1000 in incubation buffer. Incubate for 1 hour at room temperature.
 - 11. Wash as in step 4.
 - 12. Add 100 μl/well of appropriate substrate solution e.g. p-nitrophenyl-phosphate (pNPP), available from Mabtech product code 3652-P10.
 - 13. Measure the optical density (405 nm for pNPP) in an ELISA reader after suitable developing time.



Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the standards ISO 9001:2015 & ISO 13485:2016.





The products are for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages there from.

Mabtech AB (Head Office)

Sweden Tel: +46 8 716 27 00 mabtech@mabtech.com Mabtech, Inc.

USA Tel: +1 513 871-4500 mabtech.usa@mabtech.com