Human Interleukin-12/-23 (p40) Monoclonal antibody FOR NEUTRALIZATION

General


Antibody:	МТ3279Н
Product code:	3450-0N-500
Quantity:	500 µg
Immunogen:	Recombinant human IL-23
Isotype:	Mouse IgG2b
Specificity:	Native and recombinant p40 subunit of human IL-12 and IL-23
Concentration:	1 mg/ml
Buffer:	Supplied in sterile filtered (0.2 μ m) PBS. Contains no preservatives. The endotoxin content in the product is below 10 EU/mg.
Purification:	Purified from <i>in vitro</i> cultures by protein G affinity chromatography.
Storage:	Store product at 4-8°C or frozen at -20°C or below. Avoid repeated freezing/thawing.
Applications:	Neutralization of human IL-12 and IL-23 bioactivity. It is recommended to establish the optimal antibody concentration for the assay system used.

Neutralization

Neutralization assay: Neutralization of human IL-12 and human IL-23 induced IFN- γ production (see graphs). To determine the neutralizing capacity and ND₅₀ on IL-12 induced IFN- γ production, mixtures of various concentrations of mAb MT3279H, human recombinant IL-12 (10 ng/ml) and recombinant IL-1 β (50 ng/ml) were added to human PBMC. For determination of the neutralizing capacity and ND₅₀ on IL-23 induced IFN- γ production, mixtures of various concentrations of mAb MT3279H, human recombinant IL-23 (10 ng/ml) and recombinant IL-1 β (50 ng/ml) were added to human PBMC. For determination of the neutralizing capacity and ND₅₀ on IL-23 induced IFN- γ production, mixtures of various concentrations of mAb MT3279H, human recombinant IL-23 (10 ng/ml) and recombinant IL-18 (50 ng/ml) were added to human PBMC. Following overnight incubation, cell supernatants were quantified in IFN- γ ELISA.

Neutralization dose 50 (ND₅₀): MT3279H neutralized the effect of IL-12 and IL-23 at a ND_{50} value of 0.2-0.5 µg/ml when 10 ng/ml of recombinant human IL-12 or IL-23 was present.

 ND_{50} is defined as the concentration of antibody that results in 50% inhibition of the biological response under the experimental conditions specified.

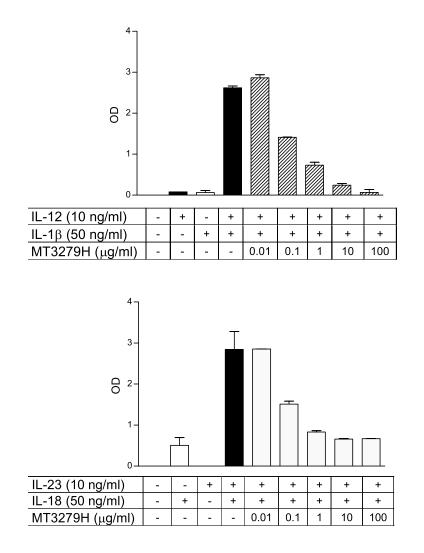



Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the standards ISO 9001:2015 & ISO 13485:2016.



The products are for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages there from.

Mabtech AB (Head Office) Sweden Tel: +46 8 716 27 00 mabtech@mabtech.com Mabtech, Inc. USA Tel: +1 513 871-4500 mabtech.usa@mabtech.com

www.mabtech.com