

# Human EBI3 ELISA development kit

Product Code: 3459-1H-6

---

CONTENTS, development kit for 6 plates:

**Vial 1 (green top)**

Monoclonal antibodies MT23A14 (300  $\mu$ l)

Concentration: 0.5 mg/ml

**Vial 2 (yellow top)**

Biotinylated monoclonal antibody MT140 (80  $\mu$ l)

Concentration: 0.5 mg/ml

**Vial 3 (white top)**

Streptavidin-Horseradish Peroxidase (80  $\mu$ l)

**Vial 4**

Recombinant human IL-27 standard

To ensure total recovery of stated quantity, vials have been overfilled.

**STORAGE:**

Shipped at ambient temperature. On arrival box 1 should be stored refrigerated at 4-8°C and box 2 should be stored frozen at -20°C.

# General


**Intended use:** For quantitative determination of native and recombinant human EBI3 in solution, e.g. cell culture supernatant. The monoclonal antibodies in this kit are specific to the EBI3 subunit. The ELISA detects both IL-27 and IL-35. The ELISA standard included is recombinant IL-27.

**Serum/plasma samples:** Please note that determination of analyte in human serum/plasma samples by this kit requires the use of ELISA diluent (product code: 3652-D2) for dilution of samples, standard and detection antibody. The diluent prevents false positive read-outs which may be caused by interference of heterophilic antibodies commonly found in human plasma and serum. The ELISA diluent has been validated using serum/plasma from normal healthy human blood donors. Please note that heterophilic antibody interference in samples from human subjects with various diseases or other conditions has not been assessed. Please contact Mabtech for further information.

**Reagents:** Antibodies are supplied in sterile-filtered (0.2  $\mu\text{m}$ ) PBS with sodium azide (0.02%). Streptavidin-HRP is supplied in PBS with 1% BSA and 0.002% Kathon CG.

**Standard range:** 15-3000 pg/ml

**Standard calibration:** No international standard exists for calibration.

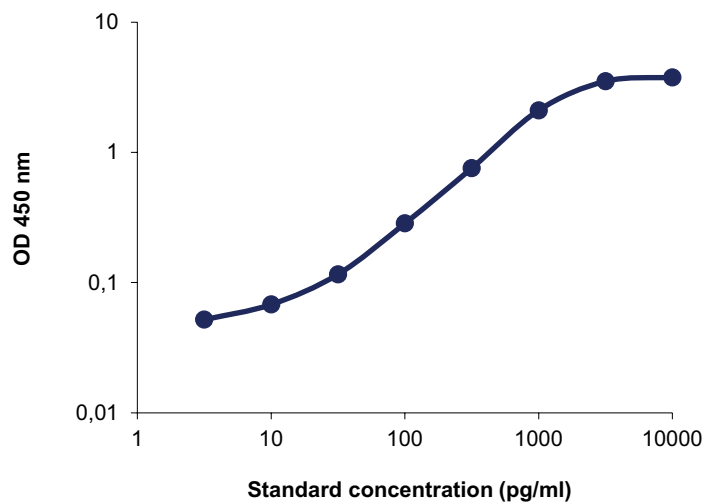

# Guidelines for Human EBI3 ELISA

- Day 1**
1. Coat a high protein binding ELISA plate with mAb MT23A14 , diluted to 2 µg/ml in PBS, pH 7.4, by adding 100 µl/well. Incubate overnight at 4-8°C.
- Day 2**
2. Wash twice with PBS (200 µl/well).
  3. Block plate by adding 200 µl/well of PBS with 0.05% Tween 20 containing 0.1% BSA (incubation buffer). Incubate for 1 hour at room temperature.
  4. Wash five times with PBS containing 0.05% Tween20
  5. Prepare standard by reconstituting contents of vial 4 in 1 ml PBS to a concentration of 0.5 µg/ml. Leave at room temperature for 5 minutes and then vortex the tube. The stock solution should be used immediately or stored in aliquots at -20°C for future use. We recommend the aliquots not be refrozen after initial use. For the test, prepare dilutions of the stock using the standard range as a guideline.
  6. Add 100 µl/well of samples or standards diluted in incubation buffer or ELISA diluent for serum/plasma samples and incubate for 2 hours at room temperature.
  7. Wash as in step 4.
  8. Add 100 µl/well of mAb MT140-biotin at 0.5 µg/ml in incubation buffer or ELISA diluent for serum/plasma samples. Incubate for 1 hour at room temperature.
  9. Wash as in step 4.
  10. Add 100 µl/well of Streptavidin-HRP diluted 1:1000 in incubation buffer. Incubate for 1 hour at room temperature. **Please note that sodium azide used in buffers will inhibit HRP activity.**
  11. Wash as in step 4.
  12. Add 100 µl/well of appropriate substrate solution e.g. TMB, available from Mabtech product code 3652-F10.
  13. Measure the optical density in an ELISA reader after suitable developing time. If required stop the reaction first.

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the standards ISO 9001:2015 & ISO 13485:2016.



**The products are for research use only.**

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages there from.

Mabtech AB (Head Office)  
Sweden  
Tel: +46 8 716 27 00  
mabtech@mabtech.com

Mabtech, Inc.  
USA  
Tel: +1 513 871-4500  
mabtech.usa@mabtech.com