

Human Apolipoprotein M ELISA^{PRO} kit

Kit for 2 plates

Product code: 3716-1HP-2

TABLE OF CONTENTS

PAGE

Kit description 4

Reagents provided 4

Storage and expiry date..... 5

Materials required but not supplied..... 5

Safety issues 5

Procedural notes 5

Preparations..... 6

Assay procedure 8

Performance of the assay..... 9

Limitations of the assay..... 10

Dilution guidelines 10

KIT DESCRIPTION

Mabtech's enzyme-linked immunosorbent assay, ELISA^{PRO}, is a complete kit for the quantification of human apolipoprotein M (apoM) in biological fluids such as serum, plasma and cell culture supernatants. The assay utilizes ELISA strip plates pre-coated with a capture monoclonal antibody (M5), to which samples are added. Captured apoM is detected by adding a biotinylated monoclonal antibody (M7-biotin) followed by streptavidin-horseradish peroxidase (SA-HRP). Addition of the enzyme substrate TMB will result in a colored substrate product with an intensity that is directly proportional to the concentration of apoM in the sample. The concentration of the apoM in the sample is determined by comparison to a serial dilution of purified apoM standard analyzed in parallel.

REAGENTS PROVIDED

Component	Quantity	Description/comments
Precoated 96-well strip plate (12 strips x 8 wells)	2 plates	Plate coated with anti-human apoM (M5). Supplied in foil bag with desiccant.
apoM standard	1 vial	Lyophilized purified apoM. See "Preparations" for reconstitution and dilution.
Detection antibody	1 vial (30 µl)	Biotinylated monoclonal antibody, anti-human apoM (M7-biotin) (0.5 mg/ml). Dilute before use.
SA-HRP	1 vial (30 µl)	Streptavidin-horseradish peroxidase conjugate. Dilute before use.
Standard reconstitution buffer AB (ready-to-use)	1 vial (1ml)	For reconstitution of lyophilized apoM standard.
Wash buffer concentrate (20x)	1 bottle (120 ml)	For all wash steps. Dilute before use.
Apo ELISA buffer concentrate 5x	1 bottle (60 ml)	Protein-containing buffer for: 1) dilution of all samples; 2) serial dilution of standard; 3) dilution of detection antibody. Dilute before use.
SA-HRP diluent (ready-to-use)	1 bottle (25 ml)	For dilution of SA-HRP.
TMB substrate (ready-to-use)	1 bottle (25 ml)	3,3',5,5' tetramethylbenzidine (TMB) enzyme substrate solution containing hydrogen peroxide.
Stop solution (ready-to-use)	1 bottle (25 ml)	0.18 M H ₂ SO ₄
Adhesive plate covers	6	To cover plates during incubations.

To ensure total recovery of stated quantity, bottles and vials are filled with larger volume than indicated.

STORAGE AND EXPIRY DATE

Shipped at ambient temperature. On arrival all components of the kit, with the exception of the lyophilized standard, should be stored at 2-8°C. **Please note that the lyophilized standard should be kept at -20°C.** After reconstitution of the lyophilized standard, it should be aliquoted and kept at -20°C. The expiry date for the unopened kit is indicated on the box. We recommend to use opened kit components within one month.

MATERIALS REQUIRED BUT NOT SUPPLIED

- Microplate reader capable of reading at 450 nm; preferably able of subtracting a reference wavelength between 570-650 nm.
- ELISA plate washer; automatic (adaptable for ELISA strip plates) or manual (e.g. multipette or squirt bottle).
- Precision pipettes and tips.
- Beakers, flasks and graduated cylinders necessary for reagent preparations.
- Tubes for standard and sample dilutions.
- Distilled or deionized water.

SAFETY ISSUES

- The stop solution, 0.18 M H₂SO₄ (<1%), is irritating to eyes and skin and should be handled with care.
- The standard should be handled with care due to unknown effects of exposure.
- Buffers and reagents in solution contain 0.002% of the preservative Kathon CG, a potential contact allergen which may cause sensitization by skin contact.
- Human and animal samples should be treated as potentially hazardous biological material.
- All material and samples should be disposed of in accordance with local regulations.

PROCEDURAL NOTES

- Do not combine components from different kit lots or components from other suppliers.
- Dilutions of standard and samples can be prepared in plastic or glass tubes.
- Total sample or standard volume added per well should not exceed 100 µl.
- The use of duplicates for each dilution of the standard, samples and blank is recommended.
- Prior to adding new reagents to the wells, ensure that there is no residual wash buffer remaining in the wells. The wells can be emptied by turning the plate upside down followed by gently tapping the plate against clean tissue paper. Please be careful to avoid that the strips fall out of the frame e.g. grip the plate by the middle.
- Serum/plasma components present in the sample may have an impact on the performance of the assay. For this reason, all samples should be diluted at least 2x in "Apo ELISA buffer".
- Sodium azide and other nucleophilic substances (often used as preservative in various buffers) interfere with the activity of horseradish peroxidase. Thus, avoid the use of other wash buffers or solutions that may contain such interfering substances.

PREPARATIONS

Plates

Allow the plates to adjust to room temperature (18-25°C) before opening the bags.

Plan the experiment to include a standard curve and an assay background control (8 x 2 wells), blank (2 wells) and sample wells. To the blank wells, only "TMB substrate" and "Stop solution" should be added. Before the analysis of absorbance values, the mean value of the blank wells should be subtracted from the standard, the assay background control and the sample values.

Assemble the required numbers of strips in the plate frame. The strips used for the experiment can be marked e.g. with a marker pen. Store the remaining strips in the foil bag containing the desiccant at 4-8°C.

Wash buffer

Prepare the required volume of wash buffer by diluting "Wash buffer concentrate" 20x with distilled or deionized water. For 1 plate, prepare 1000 ml wash buffer by adding 50 ml "Wash buffer concentrate" to 950 ml distilled or deionized water.

Apo ELISA buffer

Prepare the required volume of Apo ELISA buffer by diluting "Apo ELISA buffer concentrate 5x" 5x with distilled or deionized water. For 1 plate, prepare 150 ml Apo ELISA buffer by adding 30 ml "Apo ELISA buffer concentrate 5x" to 120 ml distilled or deionized water.

ELISA standard

Reconstitute the apoM standard in 260 μ l of standard reconstitution buffer, do not stir. It is important to wait 20 minutes before resuspending the liquid. This gives a stock solution of 0.5 μ g/ml which should be used immediately or stored in aliquots at -20°C for future use. Mix thoroughly and aliquot. Store at -20°C. Avoid repeated freeze-thaw cycles of the standard aliquotes.

Preparation of standard curve

The standard curve can be made from standard stock solution just reconstituted or from thawed aliquotes. Prepare a serial dilution of the standard no more than 30 min prior to the experiment. Duplicate wells for the standard are recommended. Dilute the standard stock solution to create a standard curve ranging from 0.027-20 ng/ml according to the scheme below. For the assay background control (0 pg/ml), use only "Apo ELISA buffer".

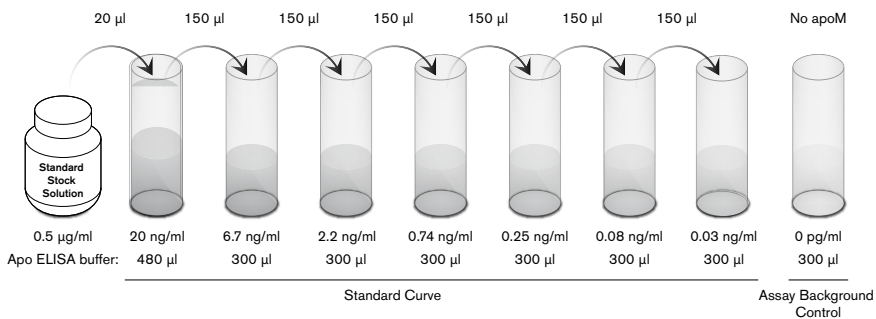


Figure 1. Recommended serial dilution of apoM standard. The volumes indicated are sufficient for duplicates. The last vial should be 0 pg/ml of standard i.e. the analyte should be omitted.

Samples

It is recommended that visible precipitate in the sample should be removed by centrifugation. Triton X-treatment of samples, necessary for apoB analysis, will not interfere with apoM analysis. Dilute all samples at least 2x in "Apo ELISA buffer" prior to use, also samples where low analyte levels are anticipated. We recommend the use of duplicate wells for all samples. It is recommended to dilute serum/plasma samples 20,000x to 40,000x. Please see dilution guidelines at page 10. Avoid repeated freezing-thawing cycles.

Detection antibody

Dilute the detection antibody "M7-biotin" 1000x in "Apo ELISA buffer" prior to use. For 1 plate, dilute 12 µl "M7-biotin" in 12 ml "Apo ELISA buffer" which will give a final "M7-biotin" concentration of 0.5 µg/ml.

SA-HRP

Dilute the "SA-HRP" 1000x in "SA-HRP diluent" prior to use. For 1 plate, dilute 12 µl "SA-HRP" in 12 ml "SA-HRP diluent".

ASSAY PROCEDURE

Throughout the assay all reagents and samples, except the "TMB substrate", should be adjusted to room temperature (18-25°C) prior to use. The "TMB substrate" should preferably be used cold.

1. Assemble the numbers of strips required for the standard curve, assay background control, the blanks and samples in the plate frame. The use of duplicates is recommended.
2. Wash the strips with 5x300 μl /well of diluted wash buffer. The wash buffer should be thoroughly removed in immediate relation to the next step.
3. Add 100 μl /well of each concentration of the diluted apoM standard and assay background control. For the samples, add 100 μl /well of pre-diluted sample. Dilution guidelines for serum/plasma samples can be found at page 10. Mix by tapping the plate.
Important! Both the serial dilution of the apoM standard and sample dilutions should be made in "Apo ELISA buffer". Leave the blank wells empty. Cover the plate with adhesive plate cover and incubate at room temperature for 2 h.
4. Wash the wells as in step 2.
5. Add 100 μl per well of the detection antibody M7-biotin, diluted 1000x to a final concentration of 0.5 $\mu\text{g}/\text{ml}$ in "Apo ELISA buffer". Leave the blank wells empty. Cover the plate with adhesive plate cover and incubate at room temp for 60 min.
6. Wash the wells as in step 2.
7. Add 100 μl /well of "SA-HRP" diluted 1000x in "SA-HRP diluent". Leave the blank wells empty. Cover the plate and incubate at room temp for 60 min.
8. Wash the wells as in step 2.
9. Develop by adding 100 μl of "TMB substrate" to all wells (including the blank wells) and incubate at room temp in the dark for 15 min.
10. Stop the color development by adding 100 μl of "Stop solution" to all wells (including the blank wells).
11. Measure the absorbance at 450 nm in a microplate reader within 15 min of the addition of the "Stop solution". If possible, use a reader capable of subtracting a reference wavelength between 570-650 nm.
12. The use of ELISA software utilizing e.g. a 4-parameter curve fitting program is recommended for the data analysis. Subtract the mean absorbance value of the blank from the standard, the assay background control and the sample values prior to creating the standard curve and determining the apoM concentrations in the samples. Note that apoM values obtained should be multiplied with the dilution factor used for each sample.

PERFORMANCE OF THE ASSAY

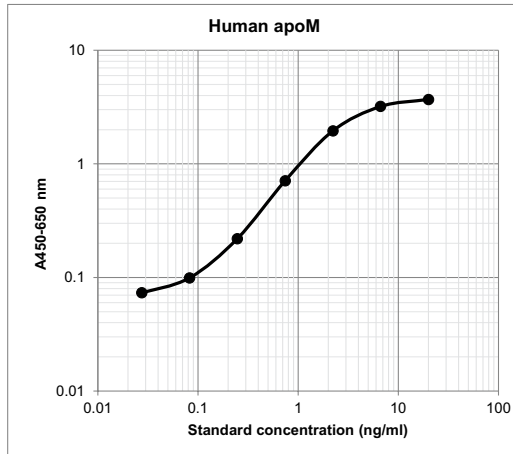


Figure 2. Display of the standard curve.

Standard range: 0.027-20 ng/ml. The standard range is the range in which determinations of analyte concentration can be done with precision, accuracy and linearity.

Sensitivity: The limit of detection of this assay is 0.01 ng/ml. It is the lowest concentration that is possible to detect but not necessarily quantify with precision and accuracy.

Dilution recovery: Dilution of serum/plasma samples gives a mid-curve recovery of 96-106% in repeated experiments.

Precision: The intraassay variation is 3.0%(CV). The interassay variation is 3.6%(CV).

LIMITATIONS OF THE ASSAY

Analysis of samples with high apoM content

The standard curve should not be extrapolated beyond the recommended standard range as these parts of the standard curve are non-linear. Samples yielding absorbance values exceeding the highest point of the standard range should be re-analyzed at a higher dilution.

Aberrant sera and plasma

The use of strongly hemolyzed and hyperlipemic serum and plasma samples may result in erroneous determinations of apoM concentrations.


Heterophilic antibodies in human serum and plasma

Heterophilic antibodies found in human serum/plasma are capable of binding to both the capture and detection antibodies used in capture ELISA. Heterophilic antibodies are found in a majority of human individuals and can, by cross-linking the assay antibodies used, result in false positive signals in capture ELISA. The "Apo ELISA buffer" provided and used for dilution of samples prevents the heterophilic antibodies from cross-linking the capture and detection monoclonal antibodies. The apoM content of serum/plasma samples can therefore be measured without interference by heterophilic antibodies. The lack of interference by heterophilic antibodies in this MABTECH kit has been validated using serum/plasma samples from normal healthy human blood donors. Please note that heterophilic antibody interference in samples from human subjects with various diseases or other conditions have not been assessed.

DILUTION GUIDELINES

Serum and plasma samples containing EDTA or heparin may be used. Samples should be aliquoted and frozen. For longer storage, -80°C is recommended. Avoid repeated freezing-thawing cycles. Before dilution, it is recommended that visible precipitate in the sample should be removed by centrifugation. Precise pipetting is important! Please change pipet tips between the different steps and ensure thorough mixing. We recommend samples to be assayed in duplicates. Sample dilutions can preferably be made in polypropylene tubes or plates. Prepare the dilutions as near as possible to the start of the assay. The indicated volumes are sufficient for duplicates. The recommended dilutions are based on serum samples from fasting healthy subjects. Samples containing levels of the analyte beyond the standard range will require other dilutions. Triton X treatment, necessary for apoB analysis, will not interfere with the analysis of other apolipoproteins.


For analysis of human ApoM

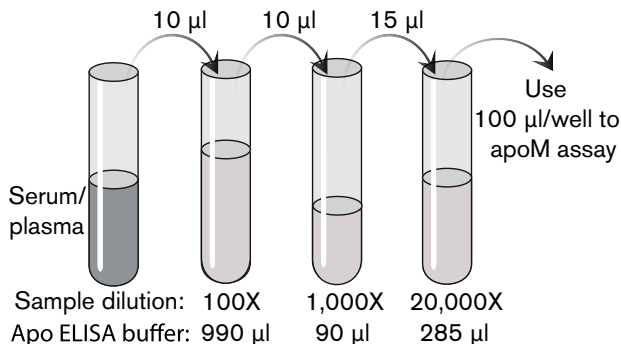

It is recommended to dilute the samples 20,000x to 40,000x. Example of dilution;

- 1) 10 µl of plasma to 990 µl Apo ELISA buffer, giving a dilution of 100x.
- 2) 10 µl 100x dilution to 90 µl Apo ELISA buffer, giving a dilution of 1,000x.
- 3) 15 µl of 1,000x dilution to 285 µl Apo ELISA buffer, giving a dilution of 20,000x.

When analyzing several apolipoproteins at the same time

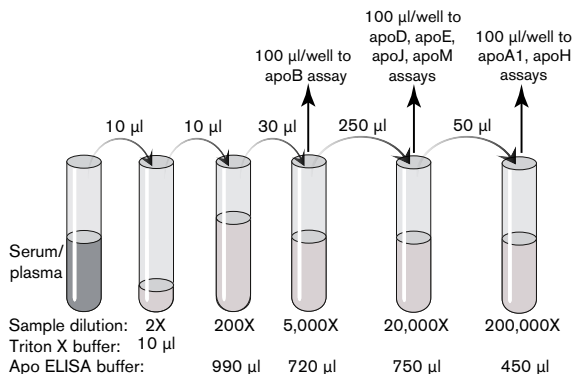

- 1) 10 µl of plasma to 10 µl of Triton X buffer, giving a dilution of 2x, vortex for 5 seconds.
- 2) 10 µl of 2x dilution to 990 µl Apo ELISA buffer, giving a dilution of 200x.
- 3) 30 µl of 200x to 720 µl Apo ELISA buffer, giving a dilution of 5,000x, use 100 µl/well to apoB ELISA
- 4) 250 µl of 5000x dilution to 750 µl Apo ELISA buffer, giving a dilution of 20,000x, use 100 µl/well to apoE, apoD, apoJ and apoM ELISAs
- 5) 50 µl of 20000x to 450 µl Apo ELISA buffer giving a dilution of 200,000x, use 100 µl/well to apoA1 and apoH ELISA

MABTECH AB
Augustendalstorget 9
Box 1233, SE-131 28 Nacka Strand
Sweden
Tel: +46 8 716 27 00
Fax: +46 8 716 27 01
E-mail: mabtech@mabtech.com

MABTECH Inc
M.E.B. 220, 3814 West Street
Cincinnati, OH 45227
USA
Toll free: +1 866 ELI-SPOT
Tel: +1 513 871 4500
Fax: +1 513 871 7353
E-mail: mabtech.usa@mabtech.com

MABTECH Australia Pty Ltd
Australia
Tel: +61 3 9470 4704
Fax: +61 3 8678 3216
E-mail: mabtech.au@mabtech.com

MABTECH AB Büro Deutschland
Germany
Tel: +49 40 4135 7935
Fax: +49 40 4135 7945
E-mail: mabtech.de@mabtech.com

MABTECH AB Bureau de liaison
France
Tel: +33 (0)4 92 38 80 70
Fax: +33 (0)4 92 38 80 71
E-mail: mabtech.fr@mabtech.com

MABTECH

NOTE; for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages therefrom.



Intertek



Intertek

Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the above standards.

For more information about Mabtech services and our products, visit

www.mabtech.com