Mouse IgA ELISA development kit

Product Code: 3835-1AD-6

CONTENTS:

Vial 1 (blue top) Monoclonal anti-IgA antibody (300 µl) Concentration: 0.5 mg/ml

Vial 2 (red top) ALP-conjugated anti-IgA antibody (150 µl)

Vial 3 Lyophilised mouse IgA standard

To ensure total recovery of stated quantity, vials have been overfilled.

STORAGE:


Shipped at ambient temperature. On arrival box 1 should be stored refrigerated at 4-8°C and box 2 should be stored frozen at -20°C.

General

Intended use: For quantitative determination of native mouse IgA in serum and plasma. May cross-react with immunoglobulins from other species.

Reagents: Monoclonal anti-IgA antibody is supplied in sterile-filtered (0.2 μ m) PBS with sodium azide (0.02%). ALP-conjugated anti-IgA antibody is supplied in 0.1 M Tris-buffer with 0.002% Kathon CG.

Recommended standard dilution: 0.1-500 ng/ml

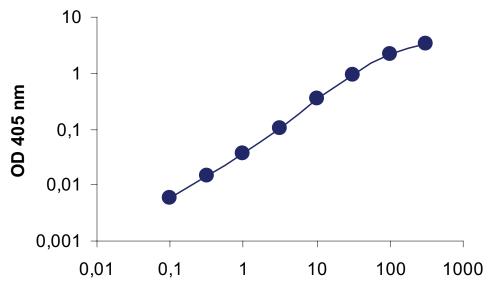

Concentration mouse IgA standard (ng/ml)

Guidelines for Mouse IgA ELISA

- **Day 1** 1. Coat a high protein binding ELISA plate with monoclonal anti-IgA antibody, diluted to 2 μg/ml in PBS, pH 7.4, by adding 100 μl/well. Incubate overnight at 4-8°C.
- **Day 2** 2. Wash twice with PBS (200 μl/well).
 - 3. Block plate by adding 200 µl/well of PBS with 0.05% Tween 20 (PBS-Tween) containing 0.1% BSA (incubation buffer*). Incubate for 1 hour at room temperature.
 - 4. Wash five times with PBS-Tween.
 - 5. Prepare mouse IgA standard by reconstituting contents of vial 3 in 500 μl PBS to make up a stock solution of 50 μg/ml. The stock solution should be used immediately or stored in aliquots at -20°C for future use. We recommend the aliquots not be refrozen after initial use. For the test, prepare dilutions of the stock using the standard range as a guideline.
 - 6. Add 100 μl/well of samples or standards diluted in incubation buffer and incubate for 2 hours at room temperature.
 - 7. Wash as in step 4.
 - 8. Add 100 μl/well of anti-IgA-ALP diluted 1:500 in incubation buffer. Incubate for 1 hour at room temperature.
 - 9. Wash as in step 4.
 - 10. Add 100 µl/well of appropriate substrate solution e.g. p-nitrophenyl-phosphate (pNPP).
 - 11. Measure the optical density (405 nm for pNPP) in an ELISA reader after suitable developing time.

* The same buffer is used for blocking and for dilution.

MABTECH AB Augustendalsvägen 19 Box 1233, SE-131 28 Nacka Strand Sweden Tel: +46 8 716 27 00 Fax: +46 8 716 27 01 E-mail: mabtech@mabtech.com

MABTECH Inc M.E.B. 220, 3814 West Street Cincinnati, OH 45227 USA Toll free: +1 866 ELI-SPOT Tel: +1 513 871 4500 Fax: +1 513 871 7353 E-mail: mabtech.usa@mabtech.com

MABTECH AUSTRALIA Pty Ltd Australia Tel: +61 3 9470 4704 Fax: +61 3 8678 3216 E-mail: mabtech.au@mabtech.com

MABTECH AB Büro Deutschland Germany Tel: +49 40 4135 7935 Fax: +49 40 4135 7945 E-mail: mabtech.de@mabtech.com

MABTECH AB Bureau de liaison France Tel: +33 (0)4 92 38 80 70 Fax:+33 (0)4 92 38 80 71 E-mail: mabtech.fr@mabtech.com

Updated on 2015-02-11

MABTECH

NOTE; for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages therefrom.



Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the above standards.

For more information about Mabtech services and our products, visit

www.mabtech.com