Human IgA ELISA development kit

Product Code: 3860-1AD-6

CONTENTS: development kit for 6 plates

Vial 1 (blue top)

Monoclonal antibody MT57 (300 $\mu l)$

Concentration: 0.5 mg/ml

Vial 2 (red top)

Monoclonal antibody MT20-ALP (80 μl)

Vial 3

Lyophilised human IgA standard

To ensure total recovery of stated quantity, vials have been overfilled.

STORAGE:

Shipped at ambient temperature. On arrival box 1 should be stored refrigerated at 4-8°C and box 2 should be stored frozen at -20°C.

General


Intended use: For quantitative determination of native human IgA in serum, plasma or cell supernatant.

Reagents: Monoclonal antibody MT57 is supplied in sterile-filtered (0.2 μ m) PBS with sodium azide (0.02%). MT20-ALP is supplied in 0.1 M Tris-buffer with 0.002% Kathon CG.

Recommended standard dilution: 0.2-100 ng/ml

Standard calibration: 1 μg of supplied standard equals 83 mU NIBSC* standard. Calibration is batch-specific.

* National Institute of Biological Standards and Control, UK.

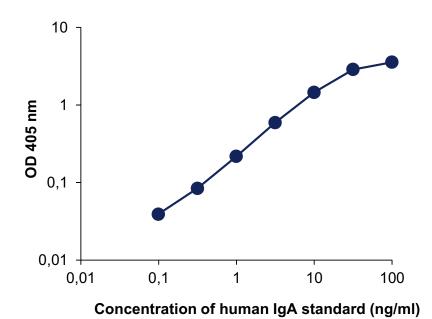

Guidelines for Human IgA ELISA

- **Day 1** 1. Coat a high protein binding ELISA plate with MT57 antibody, diluted to 2 μg/ml in PBS, pH 7.4, by adding 100 μl/well. Incubate overnight at 4-8°C.
- **Day 2** 2. Wash twice with PBS (200 μl/well).
 - 3. Block plate by adding 200 μ l/well of PBS with 0.05% Tween 20 (PBS-Tween) containing 0.1% BSA (incubation buffer*). Incubate for 1 hour at room temperature
 - 4. Wash five times with PBS-Tween.
 - 5. Prepare human IgA standard by reconstituting contents of vial 3 in 500 µl PBS to make up a stock solution of 50 µg/ml. The stock solution should be used immediately or stored in aliquots at -20°C for future use. We recommend the aliquots not be refrozen after initial use. For the test, prepare dilutions of the stock using the standard range as a guideline.
 - 6. Add 100 µl/well of samples or standards diluted in incubation buffer. It is recommended to make serial dilutions of samples to obtain results within the standard range. Incubate for 2 hours at room temperature.
 - 7. Wash as in step 4.
 - 8. Add 100 µl/well of MT20-ALP diluted 1:1000 in incubation buffer. Incubate for 1 hour at room temperature.
 - 9. Wash as in step 4.
 - 10. Add 100 μl/well of appropriate substrate solution e.g. p-nitrophenyl-phosphate (pNPP).
 - 11. Measure the optical density (405 nm for pNPP) in an ELISA reader after suitable developing time.
 - * The same buffer is used for blocking and for dilution.



Developed and manufactured by MABTECH AB, Sweden, whose quality management system complies with the standards ISO 9001:2015 & ISO 13485:2016.





The products are for research use only.

MABTECH shall not be liable for the use or handling of the product or for consequential, special, indirect or incidental damages there from.

Mabtech AB (Head Office)

Sweden Tel: +46 8 716 27 00 mabtech@mabtech.com Mabtech, Inc.

USA Tel: +1 513 871-4500 mabtech.usa@mabtech.com